

# TOTOKU

2 Megapixel  
20.1" Color Display

# CCL206

CCL Series

## 2MP color LCD with color calibration function

Dual analog/digital interfaces with a color calibration function.  
Its built-in luminance stabilizing system provides stability in luminance.



*Flat Display System  
for Medical Imaging*

**2** Megapixel  
Color Display

20.1"	Analog & Digital	High luminance 300cd/m <sup>2</sup>	Contrast 800:1	Calibration function	10Bit LUT	Multi function Stand
-------	------------------	-------------------------------------	----------------	----------------------	-----------	----------------------

Features :

▶ Color calibration function to reproduce colors faithful to the original and grayscale compliant with DICOM Part 14

By calibrating color temperature, luminance, and gamma characteristics, reproduction of colors faithful to the original and smooth display of grayscale are possible. (A calibration kit is required.)

▶ Built-in 10bit LUT (Look-up Table) to reproduce beautiful grayscale

The built-in LUT selects the most appropriate 256 shades of gray (8bit) out of 1021 shades (10bit) to realize finer and smoother grayscale.

▶ Luminance stabilizing system λ-Sentinel for stable display of medical images

Automatically corrects and luminance fluctuation due to startup drift, temperature change, and aging to provide stable luminance.

▶ Dual digital / analog interface

Analog RGB and digital DVI inputs are supported.

▶ Worldwide medical safety standards


This display is certified under various medical safety standards that are much more demanding than those for general-purpose IT equipment.

▶ Medivisor Series software (Optional)

The Medivisor Series software offers collective support for display quality control and delivers high-quality and stable medical images while realizing reduction in display quality management costs.

Specifications

Model Name	CCL206	
LCD Panel	Technology	20.1-inch, TFT Color Active matrix
	Display Area	408mm X 306mm
	Pixel Pitch	0.255mm X 0.255mm
	Contrast Ratio	800 :1 (Typ.)
	Maximum Luminance	300cd/m <sup>2</sup> (Typ.), (calibrated to 120cd/m <sup>2</sup> by factory default)
Visual Performance	Viewing Angle	178° vertical and horizontal
	Available Resolution	640 X 480 (VGA), ~1600 X 1200 (UXGA)
	Grayscale	10bit(1021) grayscale to 8bit(256) grayscale
	Color display	Full color (16.77million colors)
Input Signal	Input Sync Signal	Digital DVI / Analog RGB
	Plug and Play	DDC2B compliant
Input Power Supply	Input	AC adapter 100V ~ 240V (±10%) 50/60Hz
	Maximum Power Consumption	70W Typical (with power management feature)
Features	Calibration Control	Luminance, Gamma, Color temp. (A calibration kit is required.)
	USB Hub	USB Rev. 2.0 compliant, Self-powered USB upstream connector (x1), USB downstream connector (x2)
Approvals	UL60601-1, CSA C22.2 No. 601.1, MDD/CE, FCC-B, VCCI-B, FDA510(k)	
Physical Characteristics	Dimensions (incl. tilt stand)	Landscape : 453 (W) X 461 / 523 (H) X 220.0 (D) mm Portrait : 353 (W) X 511 / 573 (H) X 220.0 (D) mm
	Weight	Net: approx. 11.4kg, Without tiltstand: approx. 7.0kg
	Tilt stand	Tilt, Swivel, Portrait / Landscape
	Mount	100mm VESA mounting
	Security Slot	Antitheft security slot (Panel and tilt stand)
	Protective Filter	Antireflection-coated protective filter
Accessories	AC adapter, Power cord (3P), DVI cable, VGA cable, Serial cable, USB cable, Operation manual	

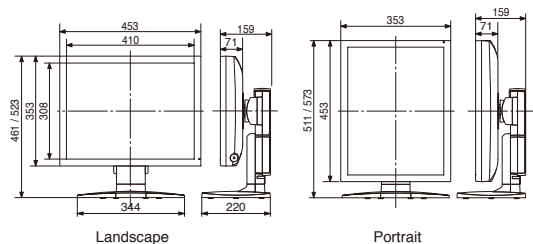

TOTOKU has obtained ISO14001 and ISO9001/ISO13485 certification which are international standards concerning environment management and quality control respectively.



**⚠ Safety Precautions** Please read the user's manual for safe and proper use.

**Do not expose the product to dust, moisture, steam, or oily smoke. It could cause fire, electric shock, or a failure.**

TOTOKU

**TOTOKU** Intelligent Devices and Solutions Dept. Sales and Marketing Division.  
TOTOKU ELECTRIC CO., LTD.  
1-3-21 Okubo, Shinjuku-ku, Tokyo 169-8543 TEL: +81 3-5273-2022 FAX: +81 3-5273-2091  
\*Osaka Branch: +81 6-6352-3541 \*Nagoya Branch: +81 52-771-1161 \*Shinetsu Branch: +81 268-34-5285

<b>USA</b> TOTOKU NORTH AMERICA, INC 1425 Greenway Dr., Ste. 565 Irving, TX 75038, U.S.A. TEL : +1-972-550-1888 E-mail : info@totoku-na.com	<b>EUROPE</b> TOTOKU Europe GmbH Hammfelddamm 6. D-41460 Neuss, Germany TEL : +49-2131-36668-0 E-mail : info@totoku.de	<b>ASIA</b> TOTOKU ELECTRIC CO., LTD. 3-21, Okubo-1 Chome Shinjuku-ku, Tokyo, Japan TEL : +81-3-5273-2053 E-mail : info-idsc@totoku.co.jp
--	---	--

Additional product information is available at <http://www.totoku.com/display/>

\*Microsoft and Windows are trademarks of the US Microsoft Corporation and are registered in the US and other countries. \*Company names and product names are the trademarks or registered trademarks of the respective companies. \*Product specifications and appearance are subject to change without notice. \*Colors in photographs may differ from actual colors due to the printing process. \*Images on screens are simulated.

Please contact the distributor below with inquiries and orders.